



Nanoindentance



Spot size



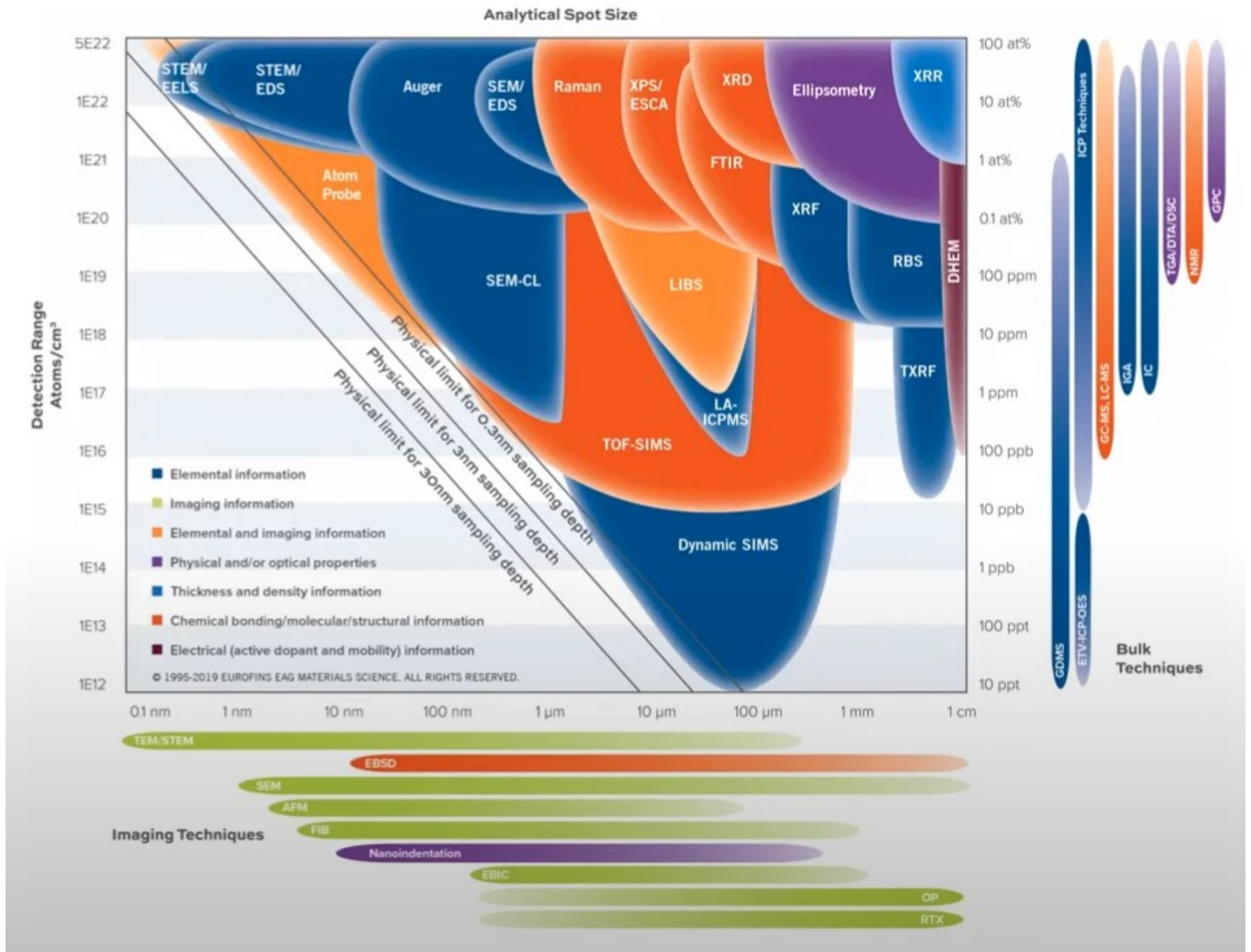


Hloubka penetrace

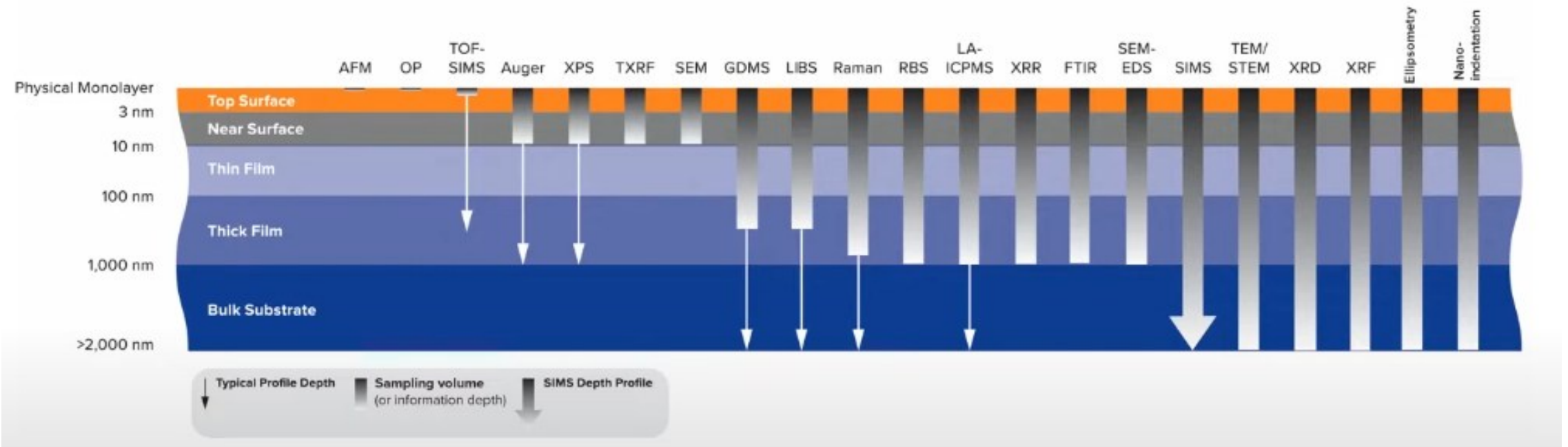


Proč nanoindentace?

Makroskopické mechanické zkoušky

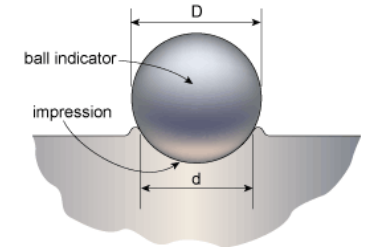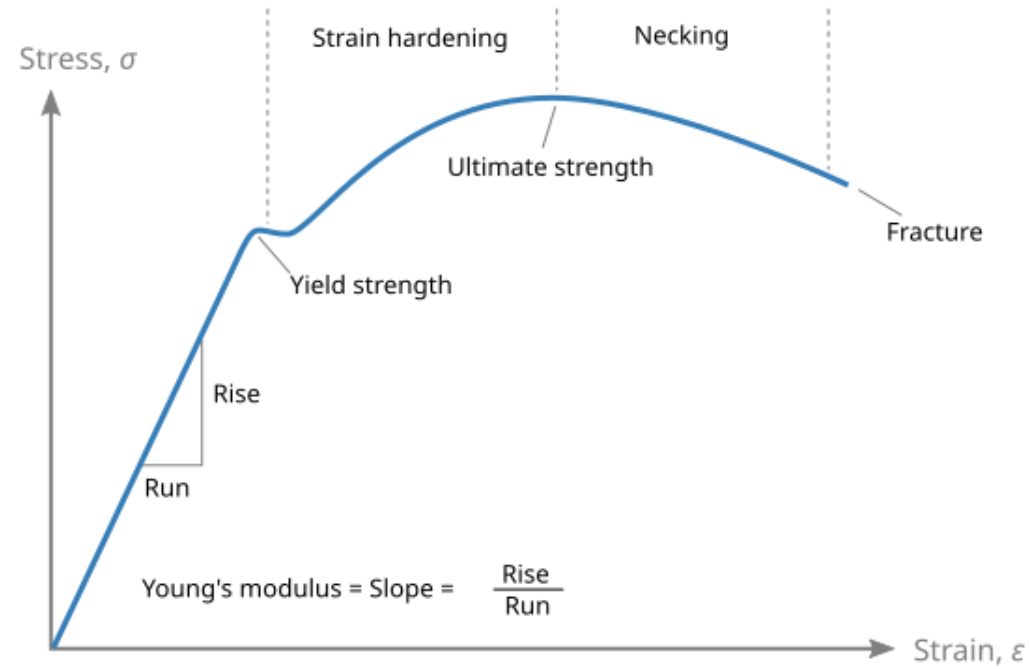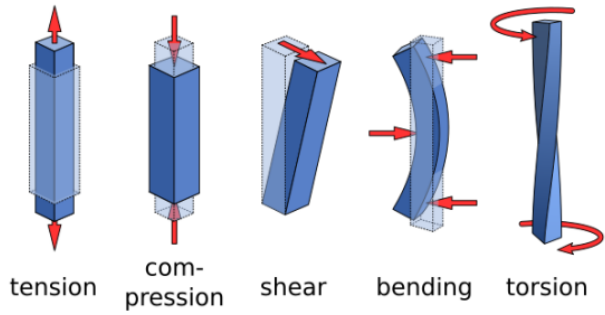


(a) Brinell indentation






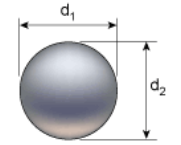


(b) measurement of impression diameter

Mikroskopické mechanické zkoušky



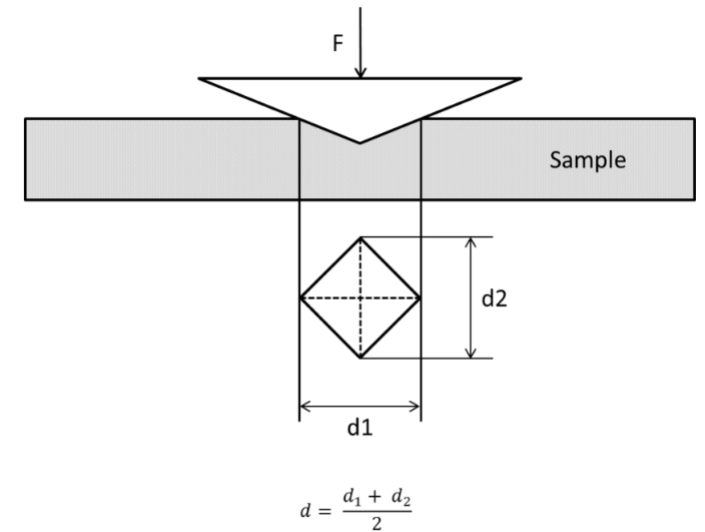


$$HV = 0.102 * \frac{2 * F * \sin \frac{136^\circ}{2}}{d^2} = 0.1891 * \frac{F}{d^2}$$

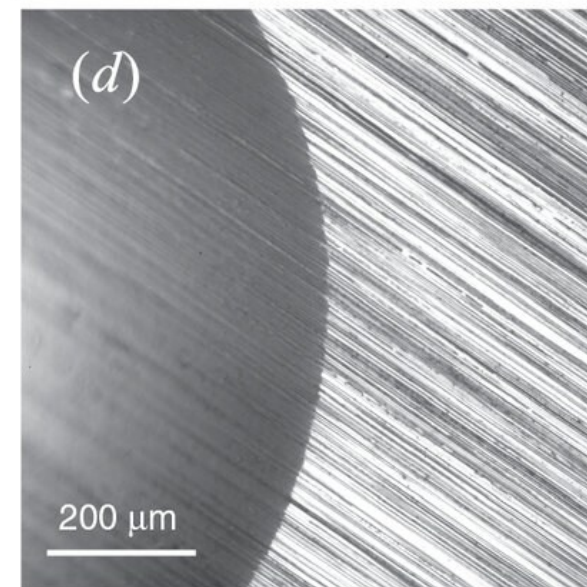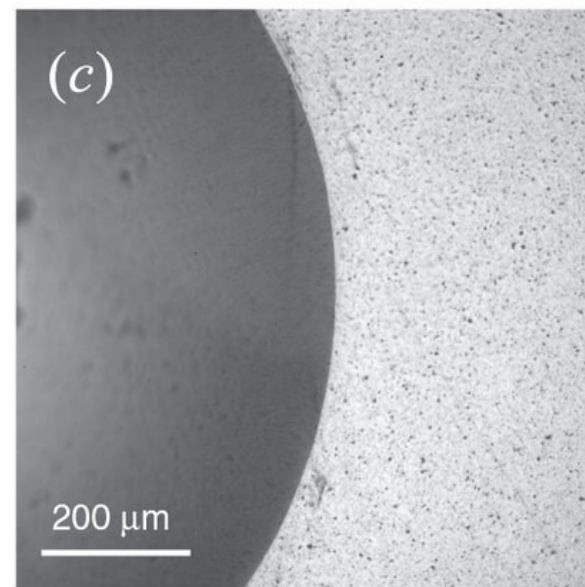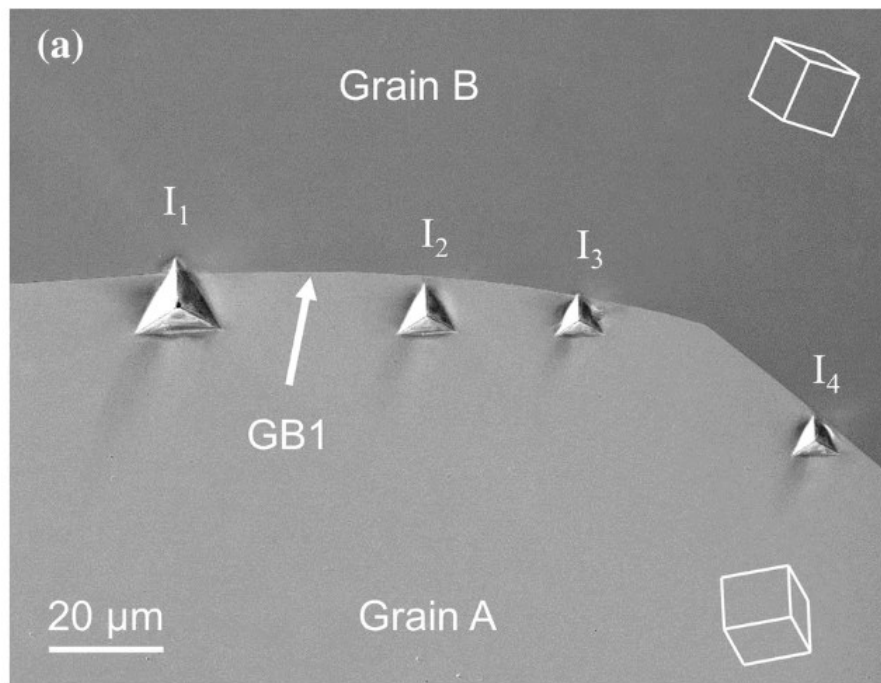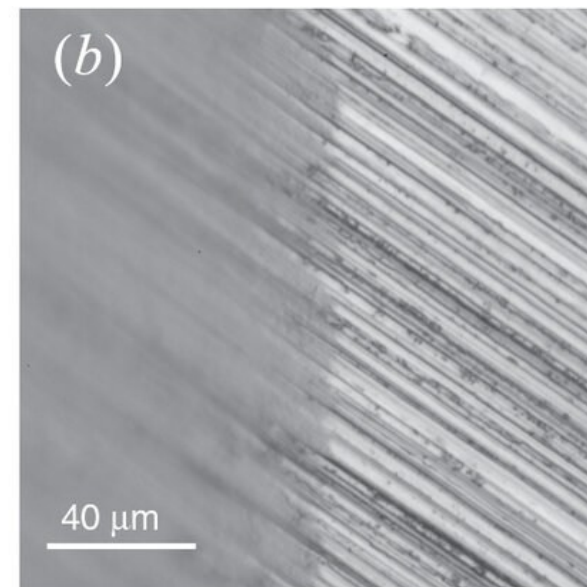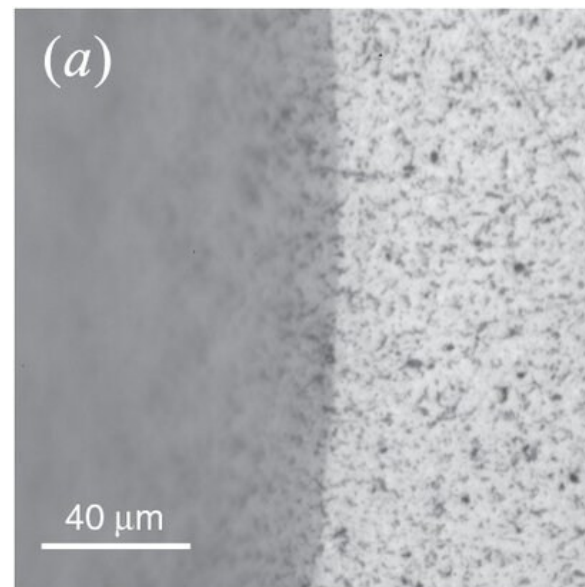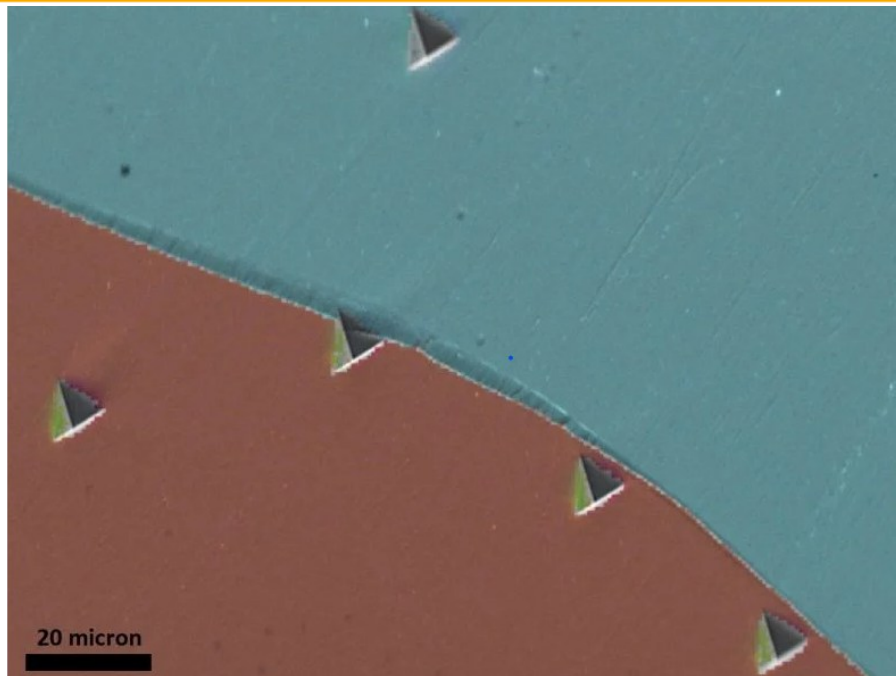


Proč nanoindentace?

Makroskopické mechanické zkoušky

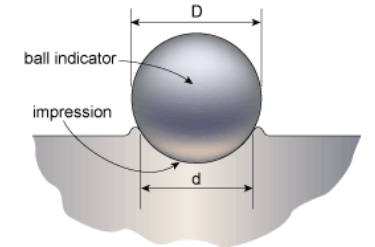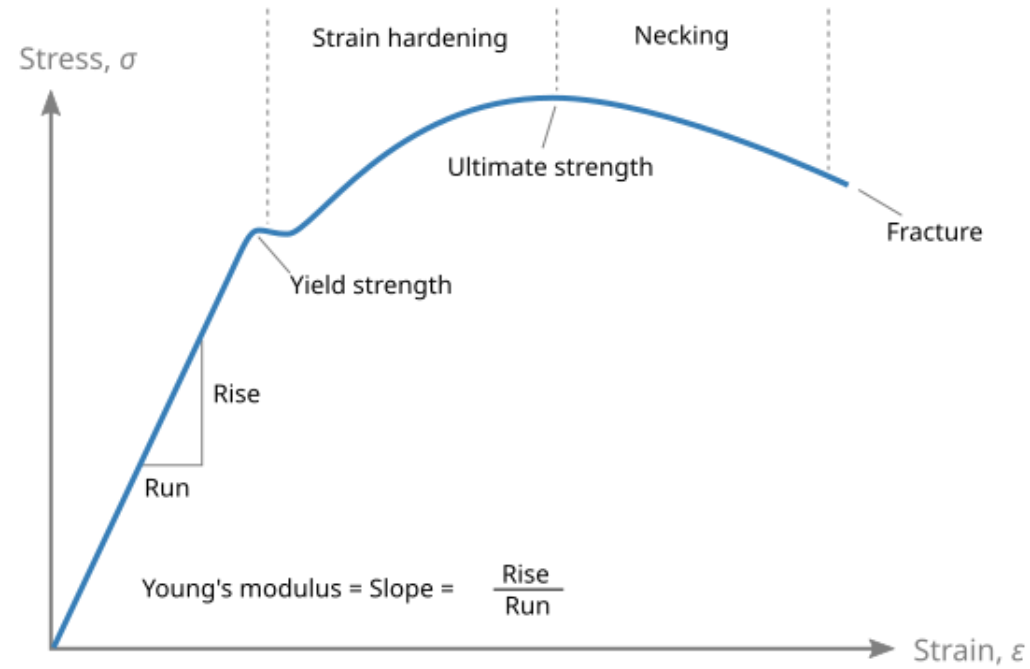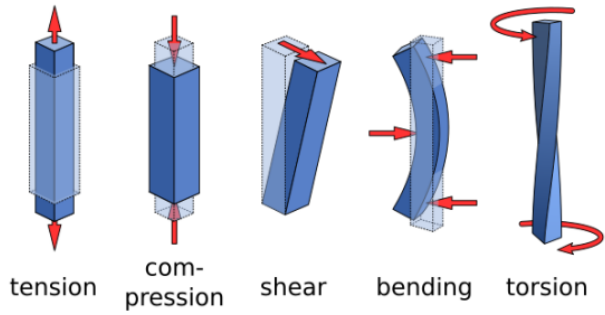


(a) Brinell indentation

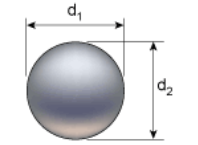

(b) measurement of impression diameter

Mikroskopické mechanické zkoušky

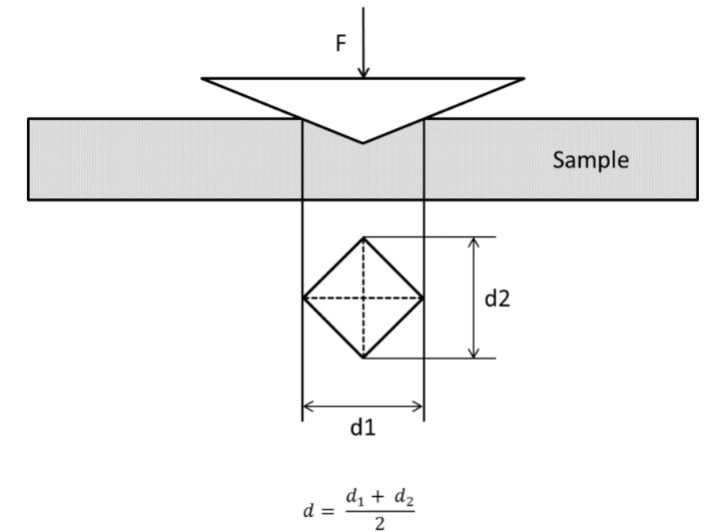

$$HV = 0.102 * \frac{2 * F * \sin \frac{136^\circ}{2}}{d^2} = 0.1891 * \frac{F}{d^2}$$


Servohydraulic Multipurpose Fatigue Test System Floor Standing with Actuator in lower base Series LFV up to 630 kN

These Tests systems meet the full spectrum of dynamic, fatigue and monotonic applications including TMF (Thermo-Mechanical Fatigue), LCF (Low-Cycle Fatigue), Fracture Mechanics, HCF (High-Cycle Fatigue), High-Strain Rate and Component Testing.

Our Servohydraulic Multipurpose LFV Series are versatile and high-performance Test Systems available in different, robust configurations covering a wide range of monotonic, dynamic and fatigue applications in the Aircraft & Aerospace, Automotive, Biomedical, Steel, Fastener, Rail or Marine industries as well Civil Engineering and R&D Testing Laboratories world-wide.

These servohydraulic systems are indispensable units of many laboratories in the field of quality control, product development or scientific research. The LFV Servohydraulic Test Systems are modular constructed, covering the full spectrum of tests requirements and can be configured to a wide range of different test setups.

Owing to over 45 years of experience in the production of a wide range of dynamic and fatigue servohydraulic testing system in different force capacities and configurations meeting the requirements of quality control, production, product research and development our machines include a numerous of features and achievements guaranteeing operational efficiency, safety and reliable testing with minimum down-time.

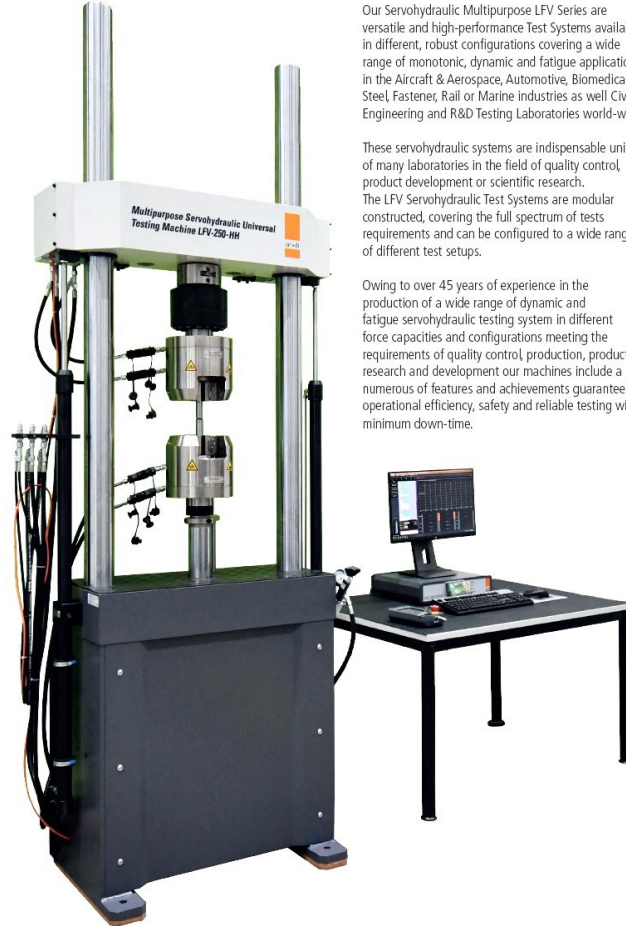

© muller-bud ag 2022. All rights reserved




Presi HZ 50-4
Vickers + Knoop



Presi Macrotech VX7

Druhy hrotů

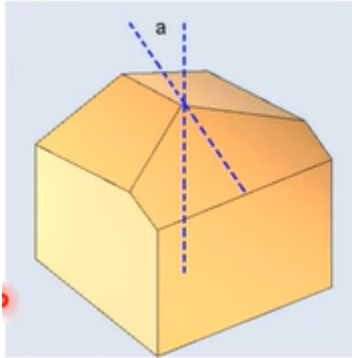


Berkovich Tip

- 3-sided pyramidal tip
- $a = 65.35^\circ$ (half angle)
- Multi-purpose tip
- Thin films >100 nm
- AIF: 150 nm radius of curvature
- **Liquid cell tip** also available


Cube Corner Tip

- 3-sided pyramidal tip
- $a = 35.26^\circ$ (half angle)
- Fine tip
- Ultra thin films & nanostructures
- Sample fracture
- Higher resolution imaging
- AIF: 40 nm radius of curvature

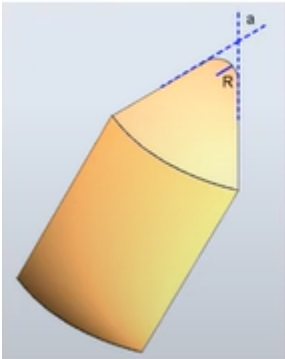

Cono-spherical Tip

- $a = 60 - 120^\circ$ (included angle)
- Very soft materials
- Lower resolution imaging
- AIF: 90° , 1 μm radius of curvature

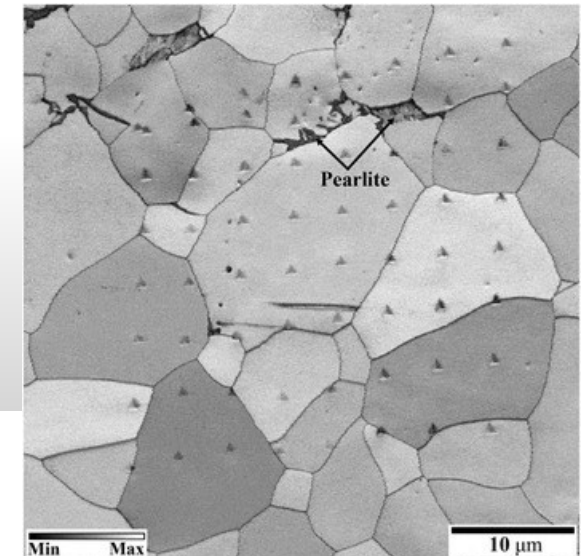

(Makroskopické zkoušky) Mechanické zkoušky



1. Materiálové zkoušky



2. Testování komponent




3. Testování stavebních materiálů



4. Kompletní a komplexní řešení




5. Speciální řešení

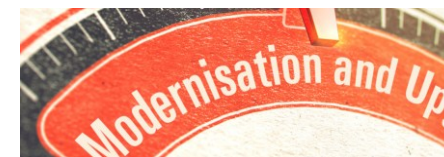

6. Modernizace, upgrade přístrojů pro mechanické testování

Mikroskopické zkoušky

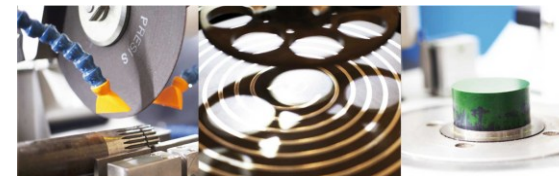

1. Krájení



2. Zalévání




3. Broušení a leštění



4. Nanoindentace




Thermo
SCIENTIFIC

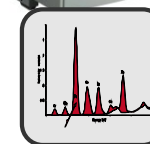

5. XRD prvková
analýza


6. Měření tvrdosti

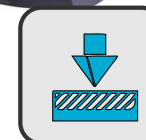

7. Odborné zaškolení
a aplikační podpora




pragolab

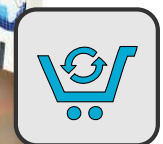

Spotřební materiál

Mikroskopické zkoušky

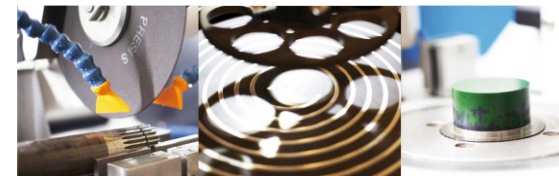

1.



pragolab

PRESI




ěštění




4. Nanoindentace

Thermo
SCIENTIFIC

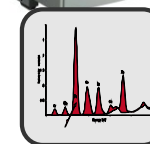

5. XRD prvková
analýza

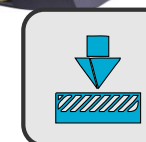

6. Měření tvrdosti




7. Odborné zaškolení
a aplikační podpora



Spotřební materiál

Mikroskopické zkoušky

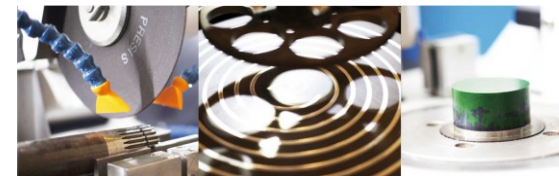


1. Krájení




4. Nanoindentace



Thermo
SCIENTIFIC

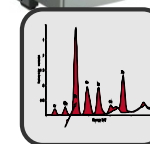


5. XRD prvková
analýza

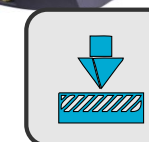

6. Měření tvrdosti




7. Odborné zaškolení
a aplikační podpora

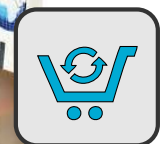

Spotřební materiál

Mikroskopické zkoušky

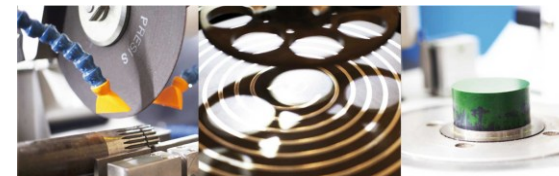


1. Krájení




2. Zalévání



Thermo
SCIENTIFIC

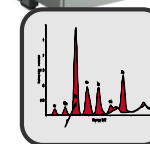

5. XRD prvková
analýza

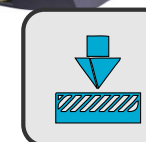


6. Měření tvrdosti




7. Odborné zaškolení
a aplikační podpora

pragolab

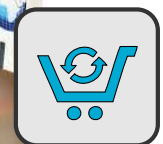

Spotřební materiál

Mikroskopické zkoušky

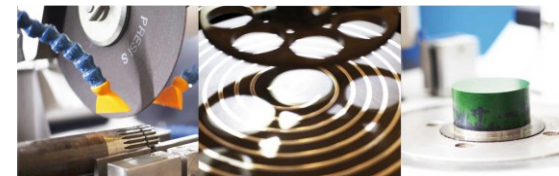

1. Krájení



2. Zalévání




3. Broušení a leštění




pragolab



Thermo
SCIENTIFIC

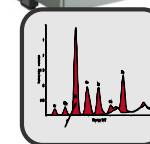

5. XRD prvková
analýza

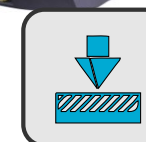

6. Měření tvrdosti




7. Odborné zaškolení
a aplikační podpora



Spotřební materiál

Mikroskopické zkoušky

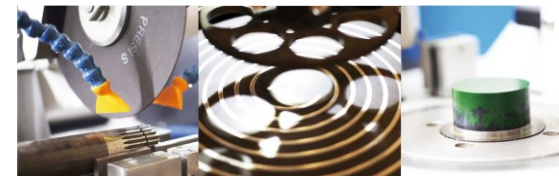

1. Krájení



2. Zalévání



3. Broušení a leštění




4. Nanoindentace




! NEW !

Thermo
SCIENTIFIC

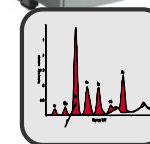

5. XRD prvková
analýza

6. Měření tvrdosti

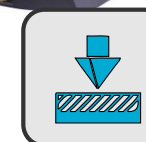

7. Odborné zaškolení
a aplikační podpora



Spotřební materiál




pragolab




Příprava vzorků

Loose objects/samples are generally not good enough to put in the nanoindenter


- **Recommended:** Super glue the sample to a magnetic AFM specimen disk (shown right)
- Super glue to a glass slide or clean wafer
- Smooth samples (wafers, glass slide) can be adhered via air suction without gluing

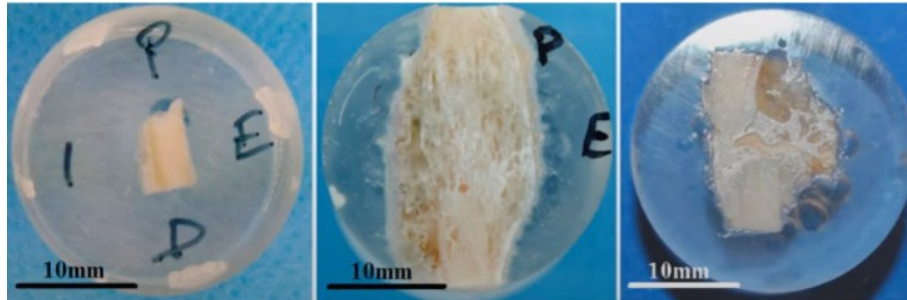

Source: Mora-Macias et. al. (2017)
<https://doi.org/10.1016/j.jmbbm.2017.05.031>
JMBBM Vol. (Oct 2017) Vol. 74, pp.236-244

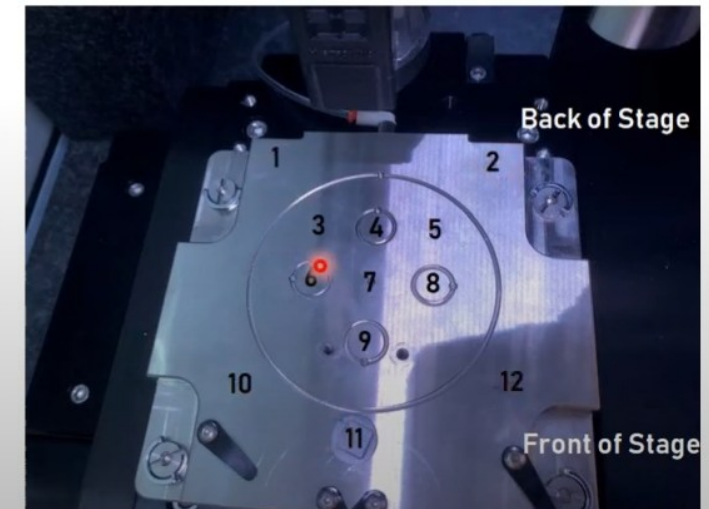

Profil povrchu a penetrační hloubka



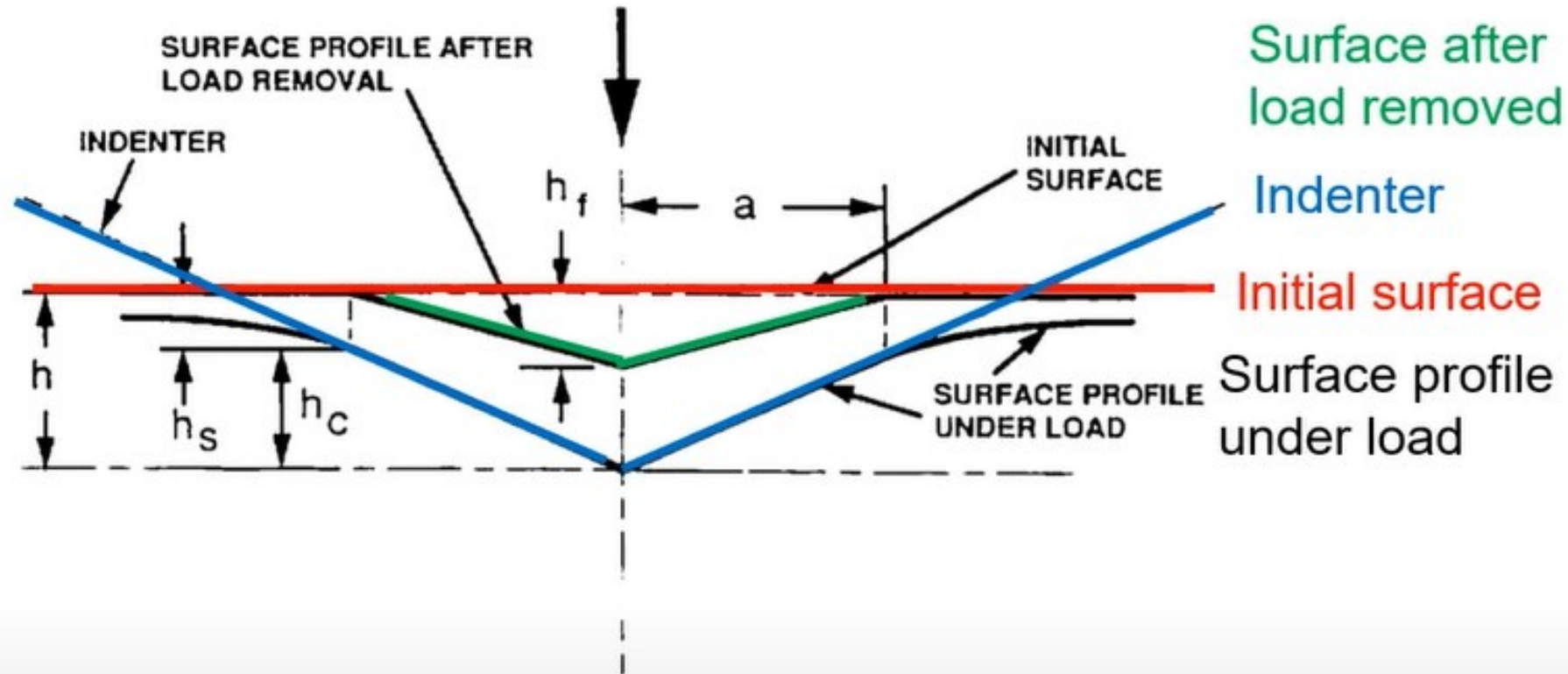


h = **penetration depth**

h_c = **contact depth** = vertical distance along which contact is made

h_s = surface displacement at perimeter of contact (**sink-in**)



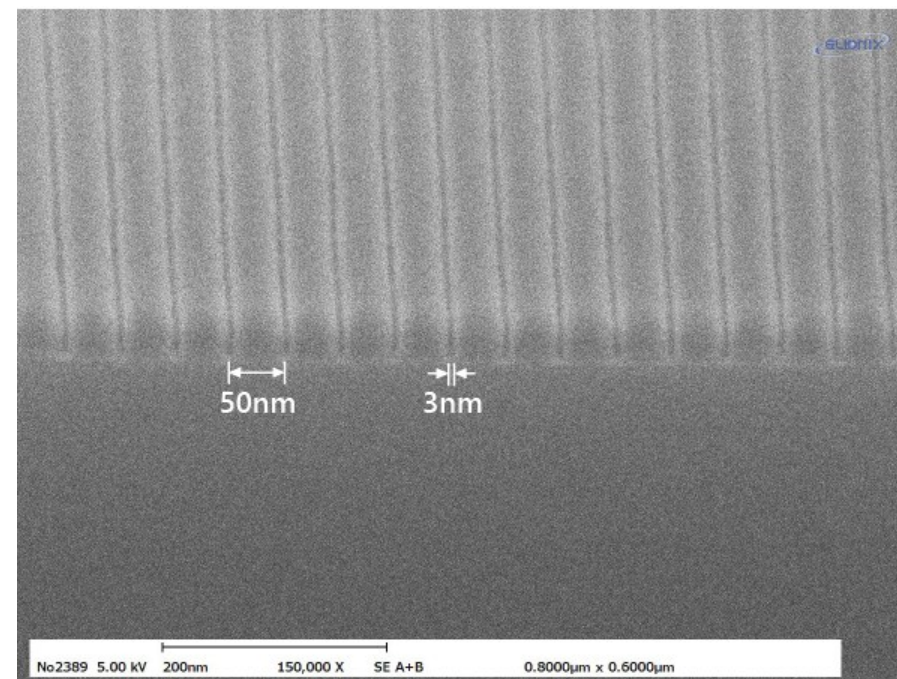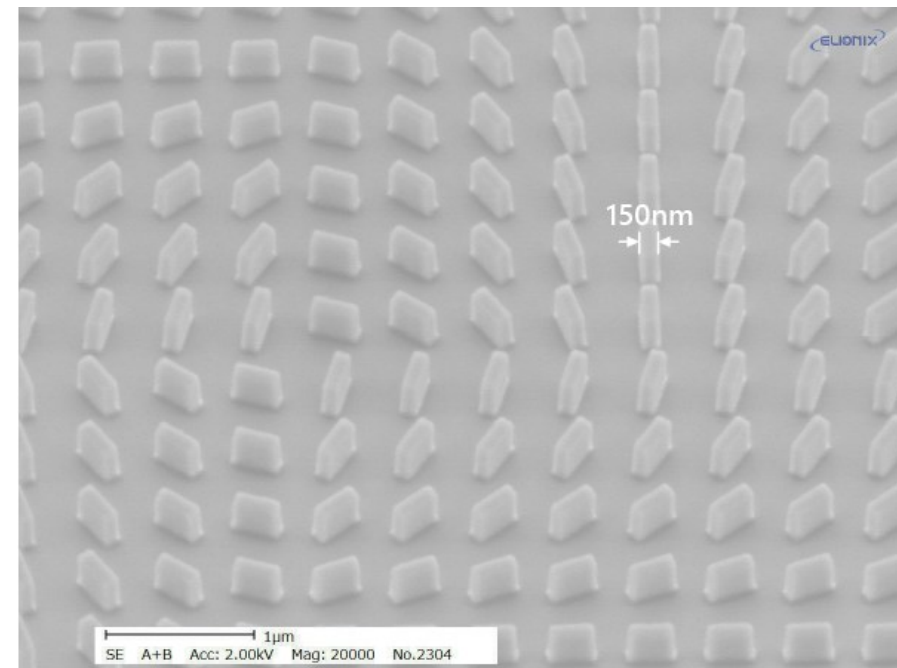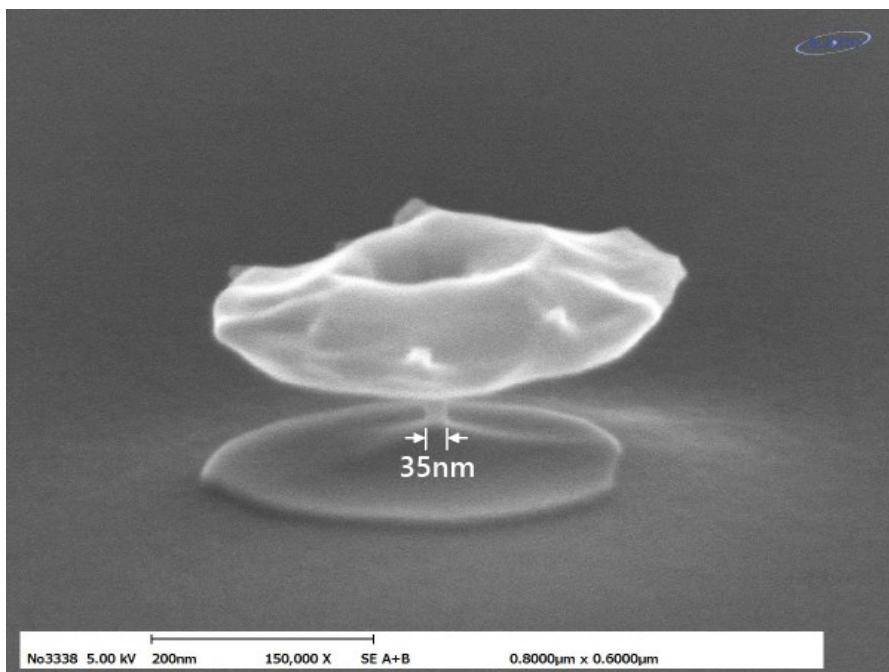



NANO INDENTATION TESTER
ENT-5

ELIONIX
Action for Innovation

Nanoindentation tester

1. Temperature control

- Prevent thermal expansion of samples and testing machines
- Temperature control in the environmental isolation shield at $\pm 0.1^{\circ}\text{C}$
- Low thermal expansion material (Nobinite) is used for machine frame and Sample holder /stage
- Reduce impact from environmental airflow

2. High precision positioning stage

- Adopts in-house developed high-precision stage technology from Elionix Electron Beam Lithography system
- 0.1 μm step in positioning
- Observe the material with magnification of 2,000x and set the measurement point.
(standard objective lens 20x and digital zoom)


3. Active vibration isolation

- Vibration-resistant wedge-shaped high-precision stage
- Active vibration isolation table



NANO INDENTATION TESTER

ENT-5

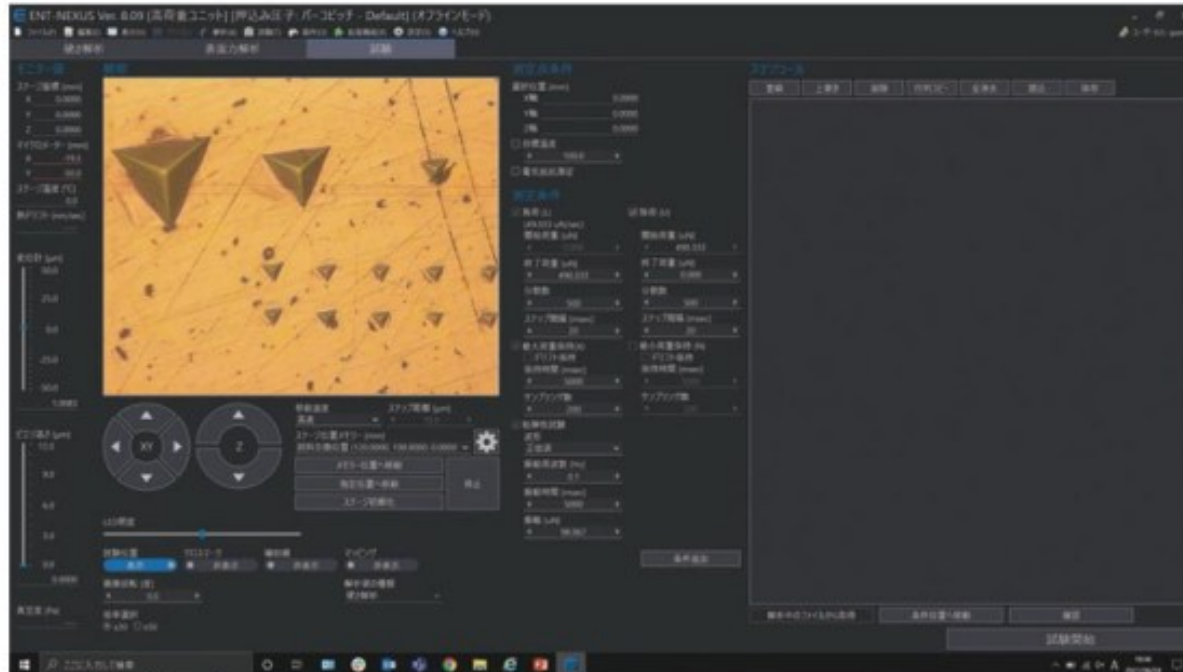

※It may differ from the actual product software

4. Software

- Pop-up alerting for daily check : indenter tip correction, temperature drift correction.
- Easy-to-use operability

5. Data reproducibility

- Improved data reproducibility by suppressing the disturbances from vibration and temperature changes
- Stable data even with continuous measurement

6. Maintainability

- Easy replacement of loading unit
- Minimum downtime with easy-change loading units


Features

● Applicable to variable kind of materials

- Plating, hard-coat thin layer
- Resin, polymer thin film
- Ultra-thin DLC coating
- Functional resin, surface modification layer
- Fine particles and powder materials

● Nano indentation test

- Mechanical characteristics such as hardness and elastic modulus on thin film or surface layer can be obtained.
- Analyze the load/unload curves to determine the characteristics. No need to observe indentations
- Compliant with ISO14577-1 / JIS Z 2255

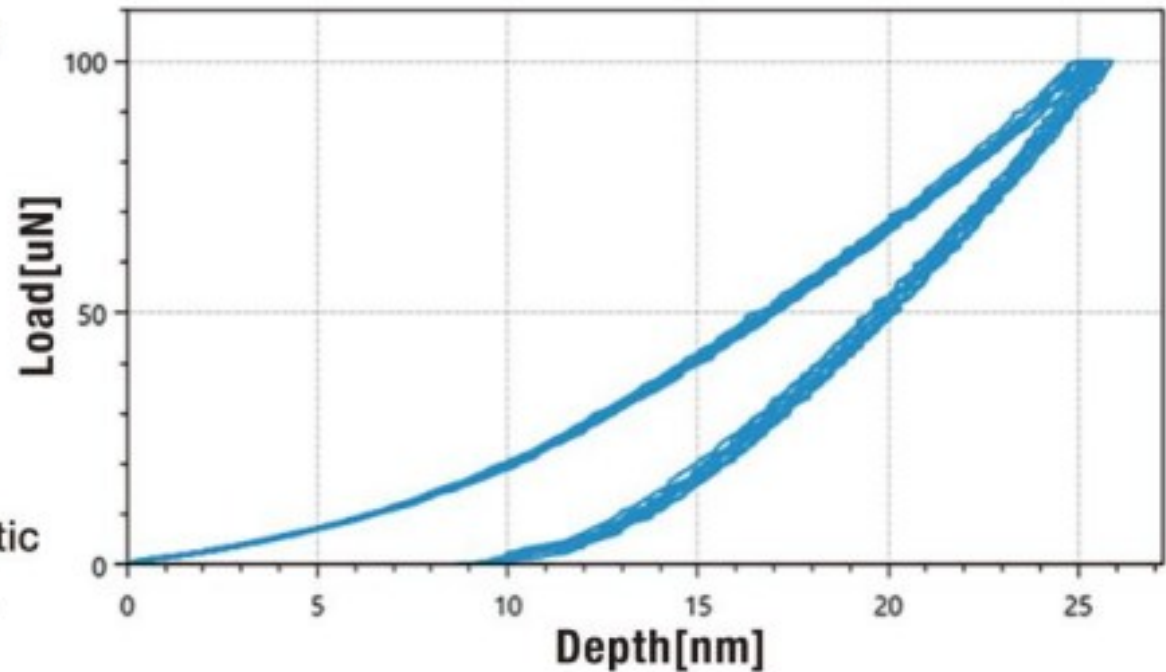

● Measurable materials



Polymer films



Electronic components
(IC, LSI, etc.)



Rubber,
Elastomer



Glass,
Ceramic



Metal



Plastic



Various fibers,
Carbon fiber



Thin film / Coating




h_c = **contact depth** = vertical distance along which contact is made

Correction for sink-in:

$$h_c = h_{max} - \varepsilon P_{max}/S$$

ϵ = geometric constant = 0.75 (paraboloid)
= 0.72 (conical)
= 1 (flat punch)

Deformace v průběhu zatížení

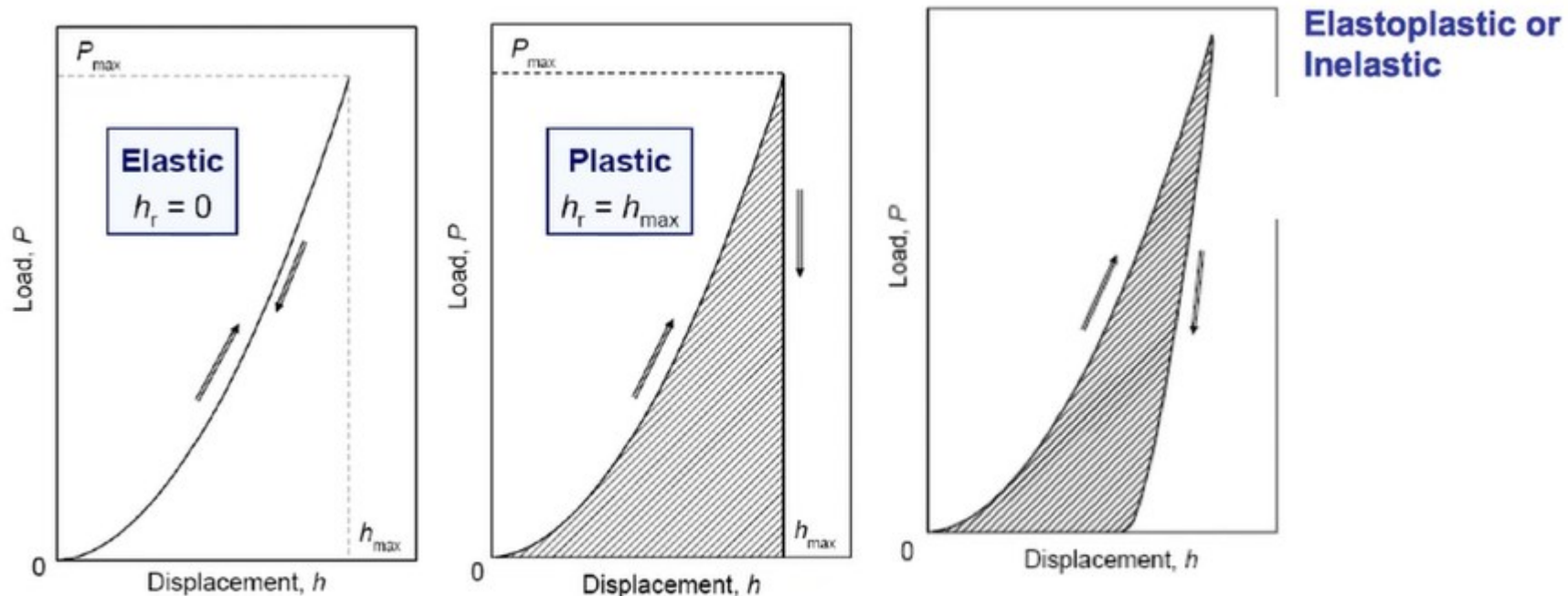


$h_r = h_f$ = residual / final depth

Measurement Example


● Loading units

High load unit ////////////////

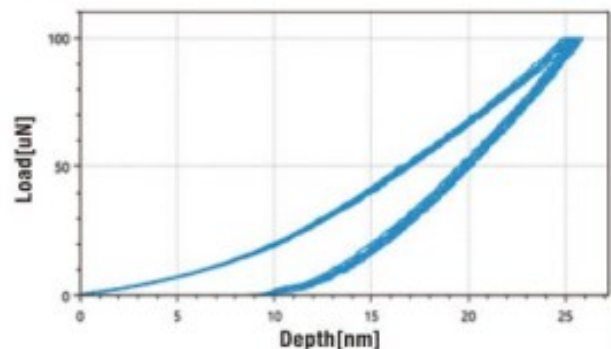

High load unit Measurement example [Fused silica]
(Load : 100 μ N Overlay 9times)

Low load unit ////////////////

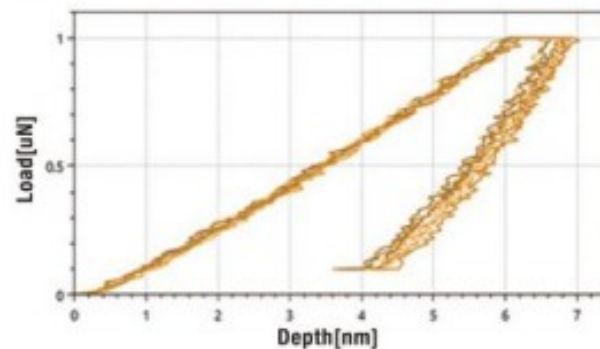

Low load unit Measurement example [resist]
(Load : 1 μ N Overlay 9times)

● Sample holder

Optional holders ////////////////



Standard sample holder (ϕ 50mm)



Heating holder
* Sample size
10 $\times$ 15 $\times$ t3mm
* Temperature range
R.T $\sim$ 250 $^{\circ}$ C (High load)
R.T $\sim$ 180 $^{\circ}$ C (Low load)



Cooling holder
* Sample size
10 $\times$ 10 $\times$ t1mm
* Temperature range
-30 $^{\circ}$ C $\sim$ 100 $^{\circ}$ C



Sample holder for Resin
embedding sample
* Sample size
Diameter 25 $\sim$ 40mm
Height 9 $\sim$ 14mm

*Customized holders to fix samples with special shapes are available.

● Measurement example

Dynamic Instrumented Indentation


- This is the method measuring continuous stiffness distribution along the depth. Applying small oscillation during loading course and calculate hardness and elastic modulus.
- Sample
Fused silica (blue), DLC coating (red)
- Test results
Fused silica shows almost constant hardness regardless of the indentation depth.
On the other hand, the DLC sample shows gradual changing in the hardness from the surface to the base material.


Underwater test (optional) ////////////////

- Hardness test in the water/solution with immersion objective lens and indenter.




Underwater sample
holder
* Sample size
 ϕ 20mm $\times$ t3mm
* Schematic

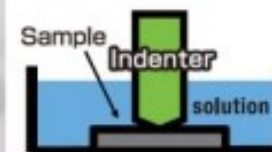

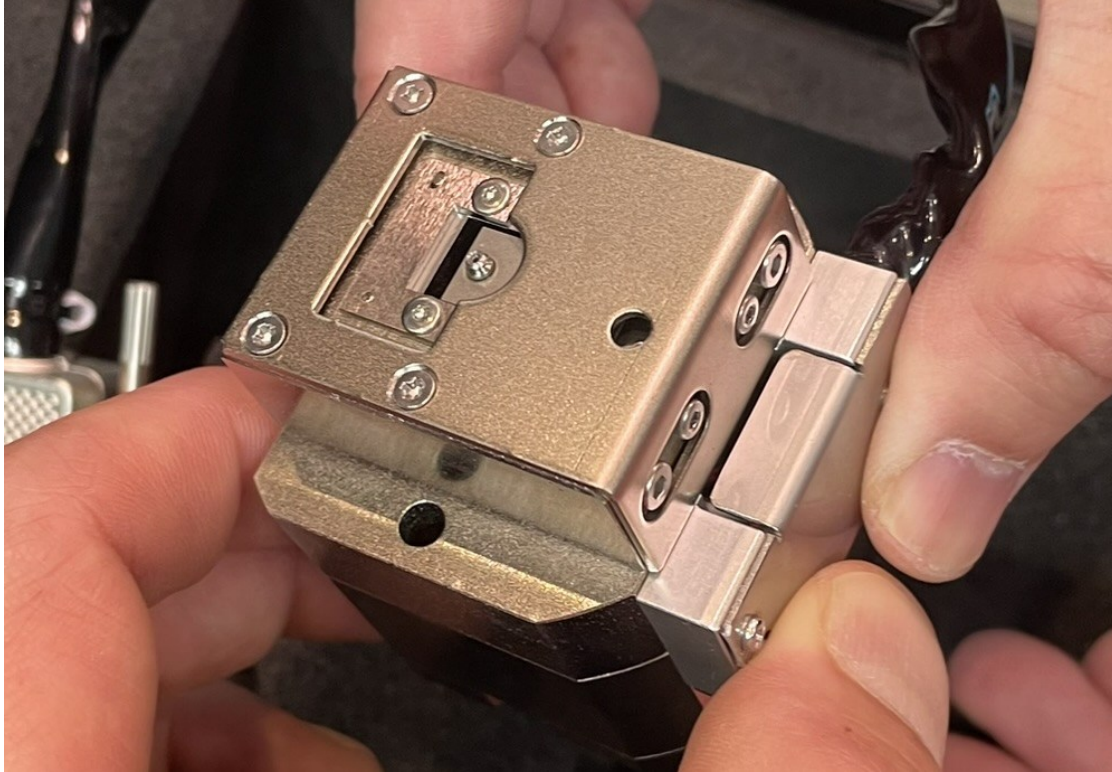


Možné hlubší poznání

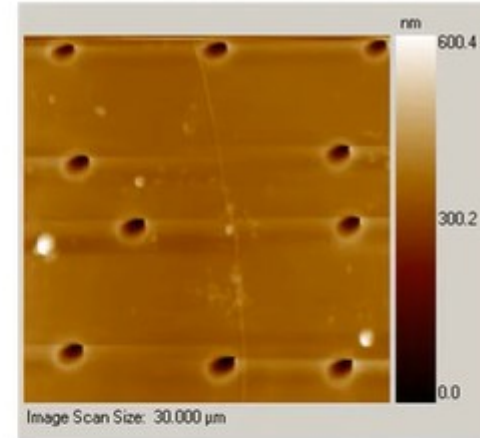

Scanned Probe Microscopy (SPM) Image of acrylic sample post-indentation


	Reduced Modulus (GPa)	Hardness (MPa)
Acrylic Untreated	5.20 GPa ± 0.03	356.5 MPa ± 6.09
Acrylic Disinfectant A	4.97 GPa ± 0.03	335.1 MPa ± 1.53
Acrylic Disinfectant B	5.28 GPa ± 0.05	355.3 MPa ± 3.26
Acrylic Disinfectant C	4.50 GPa ± 0.04	256.0 MPa ± 3.26

Average modulus and hardness measured for each sample

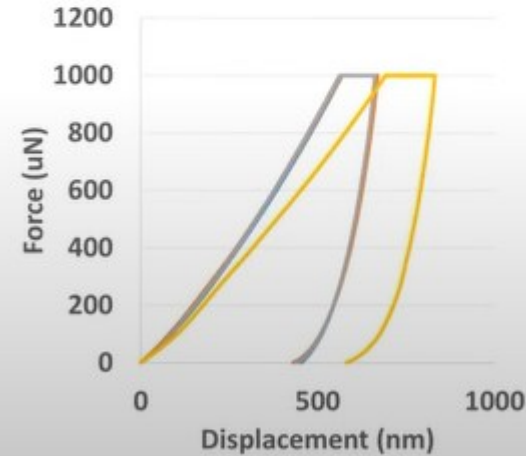


- Acrylic, Untreated
- Acrylic, soaked 24h in Disinfectant A
- Acrylic soaked 24h in Disinfectant B
- Acrylic, soaked 24h in Disinfectant C

Single force-displacement curve of each sample overlaid as measured in this study

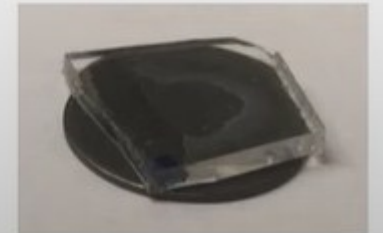


Source:
Strader et al.. **Nano Life**; Vol. 9, Iss. 4, (Dec 2019).

Test (standard function) //

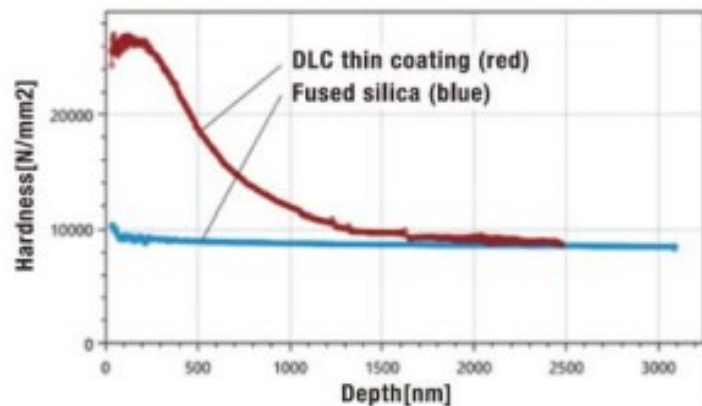

Heating / cooling viscoelasticity test (optional) //

- What is a viscoelasticity test?
This is the method to obtain elastic modulus (E'), the loss elastic modulus (E''), and the loss coefficient ($\tan\delta$). Applying small oscillation at the maximum loading.
- Sample
Ultra High Molecular Weight Poly Ethylene (PE-UHMW)
- Test results
The result with low-load unit shows same tendency as tensile test.

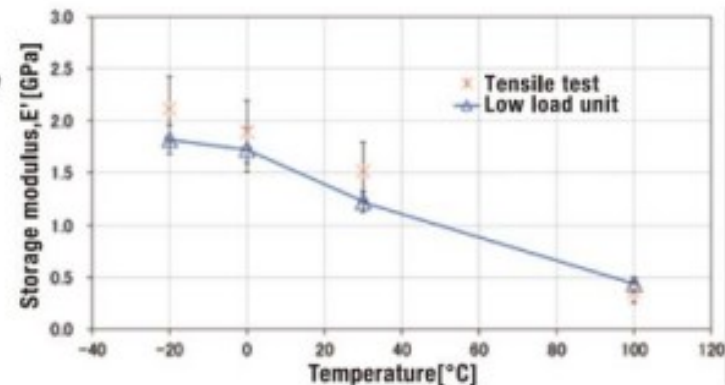


Compressive strength of micro particles test (optional) //




- This is the method to measure the fracture strength /deformation strength of particles in which a load is continuously applied by a flat punch.

JIS Z 8844:2019 Test method of fracture and deformation strength of a fine particle


* Schematic

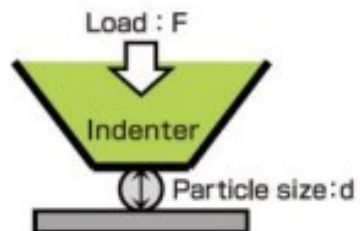

- Sample
Glass (Particle size $10\mu\text{m} \pm 0.5\mu\text{m}$)
Polystyrene (Particle size $10\mu\text{m} \pm 0.5\mu\text{m}$)
Copper (Particle size $10\mu\text{m} \pm 2\mu\text{m}$)
- Test results
In ductile polystyrene and copper, the displacement changes continuously with the load.
In brittle glass, discontinuous curve is seen when the glass is broken.

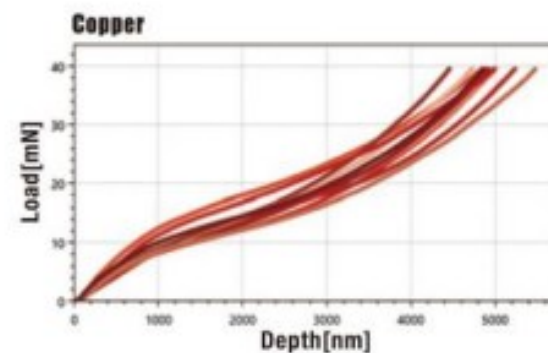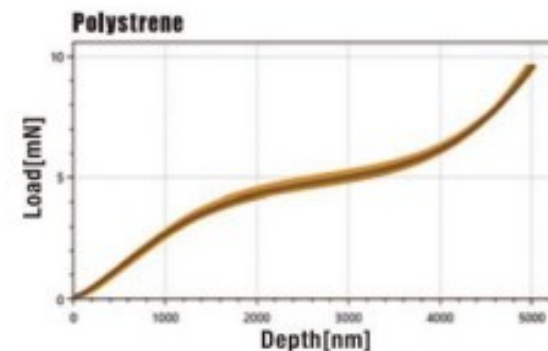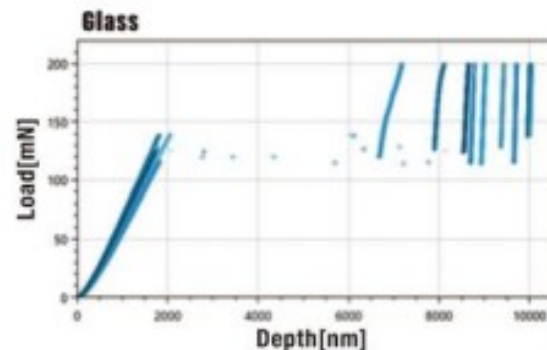

Video

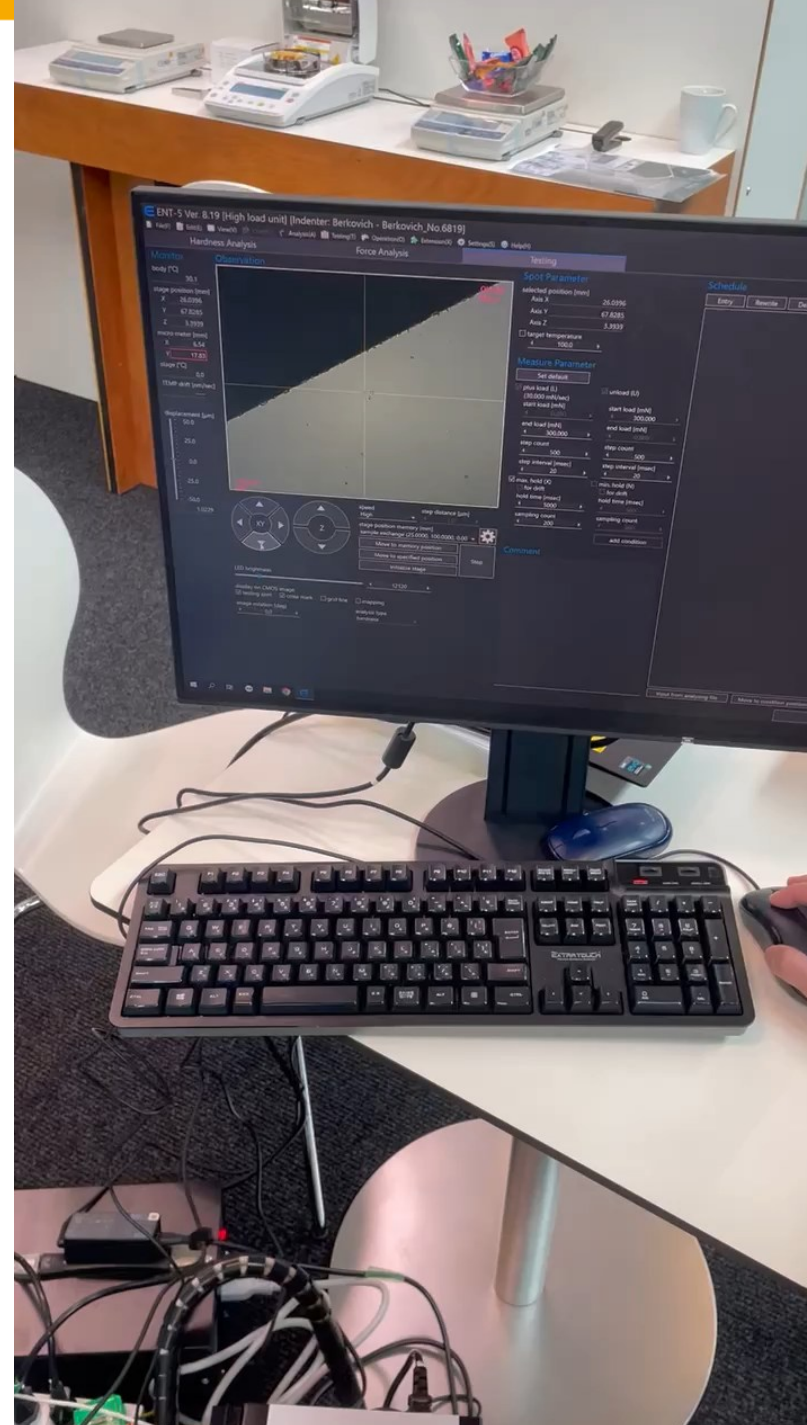

Analysis

Testing

Monitor

stage position [mm]

X 195.9152

Y 87.7779

Z 5.1695

stage [°C]

0.0

TEMP drift [nm/sec]

0.9933

displacement [μm]

50.0

25.0


0.0

-25.0

-50.0

0.9933

Observation

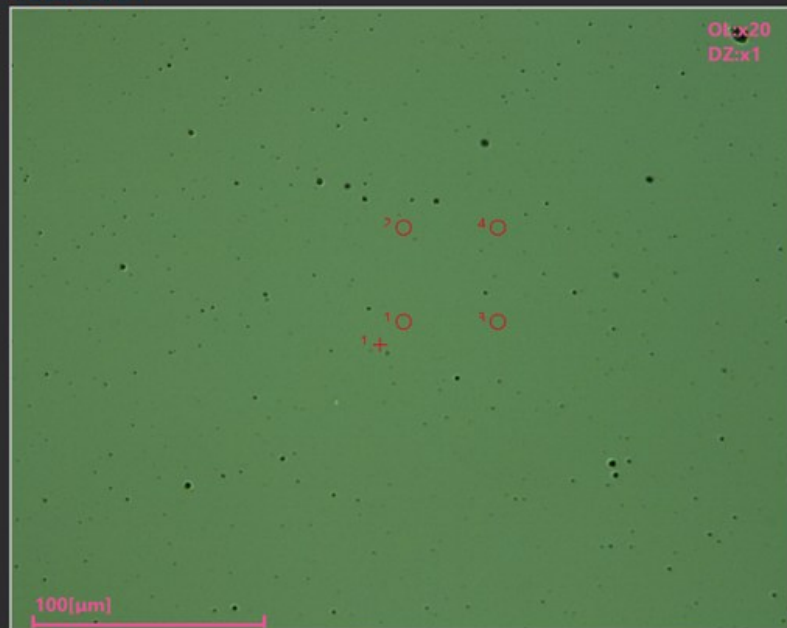

speed User step step distance [μm]

stage position memory [mm]

sample exchange (120.0000, 100.0000, 0.0)

Move to memory position

Move to specified position

Initialize stage

Stop

LED brightness

testing spot

show

cross mark

hide

grid line

hide

mapping

hide

image rotation [deg]

0.0

analysis type

硬さ解析

Spot Parameter

selected position [mm]

Axis X 195.9152

Axis Y 87.7779

Axis Z 5.1695

target temperature

100.0

Measure Parameter

plus load (L)

(1.000 mN/sec)

start load [mN]

0.000

end load [mN]

10.000

step count

500

step interval [msec]

20

max hold (X)

for drift

hold time [sec]

5.000

sampling count

200

unload (U)

start load [mN]

10.000

end load [mN]

0.000

step count

500

step interval [msec]

20

min hold (N)

for drift

hold time [sec]

5.000

sampling count

200

add condition

Comment

Schedule

Entry

Rewrite

Delete

Matrix

All clear

Load

Save

- spot001 (195.9152, 87.7779, 5.1695)[mm]
measure001 load-unload (load(0.00-10.00[mN]), max hold(5[sec]), unload(10.00-0.00[mN]))
- spot002 (195.9152, 87.8079, 5.1695)[mm]
measure001 load-unload (load(0.00-10.00[mN]), max hold(5[sec]), unload(10.00-0.00[mN]))
- spot003 (195.9452, 87.7779, 5.1695)[mm]
measure001 load-unload (load(0.00-10.00[mN]), max hold(5[sec]), unload(10.00-0.00[mN]))
- spot004 (195.9452, 87.8079, 5.1695)[mm]
measure001 load-unload (load(0.00-10.00[mN]), max hold(5[sec]), unload(10.00-0.00[mN]))

Account Setting

Basic Settings

language

English

color theme

Dark

default image format for saving

Png

character encoding form for csv file

UTF-8

new password

new password (confirm)

data directory

D:\ELIIONIX\ENT8\userdata\guest

Set default

Load guest user's parameters

Load another user's parameters

Email Notification Settings

Notify by email when exam is completed

Mail sending server Name(SMTP) Mail sending server Port Number(SMTP)

Authentication user name of mail sending server Authentication user password of mail sending server

Sender's email address

Sender's display name

Destination email address

Keywords to add to subject

Send test

Apply

Cancel

Input from analyzing file


Move to condition position

Confirm


Start testing

REPORT OUTPUT (Original data)


REFERENCE only.
Do NOT Distribute,
just to show only.

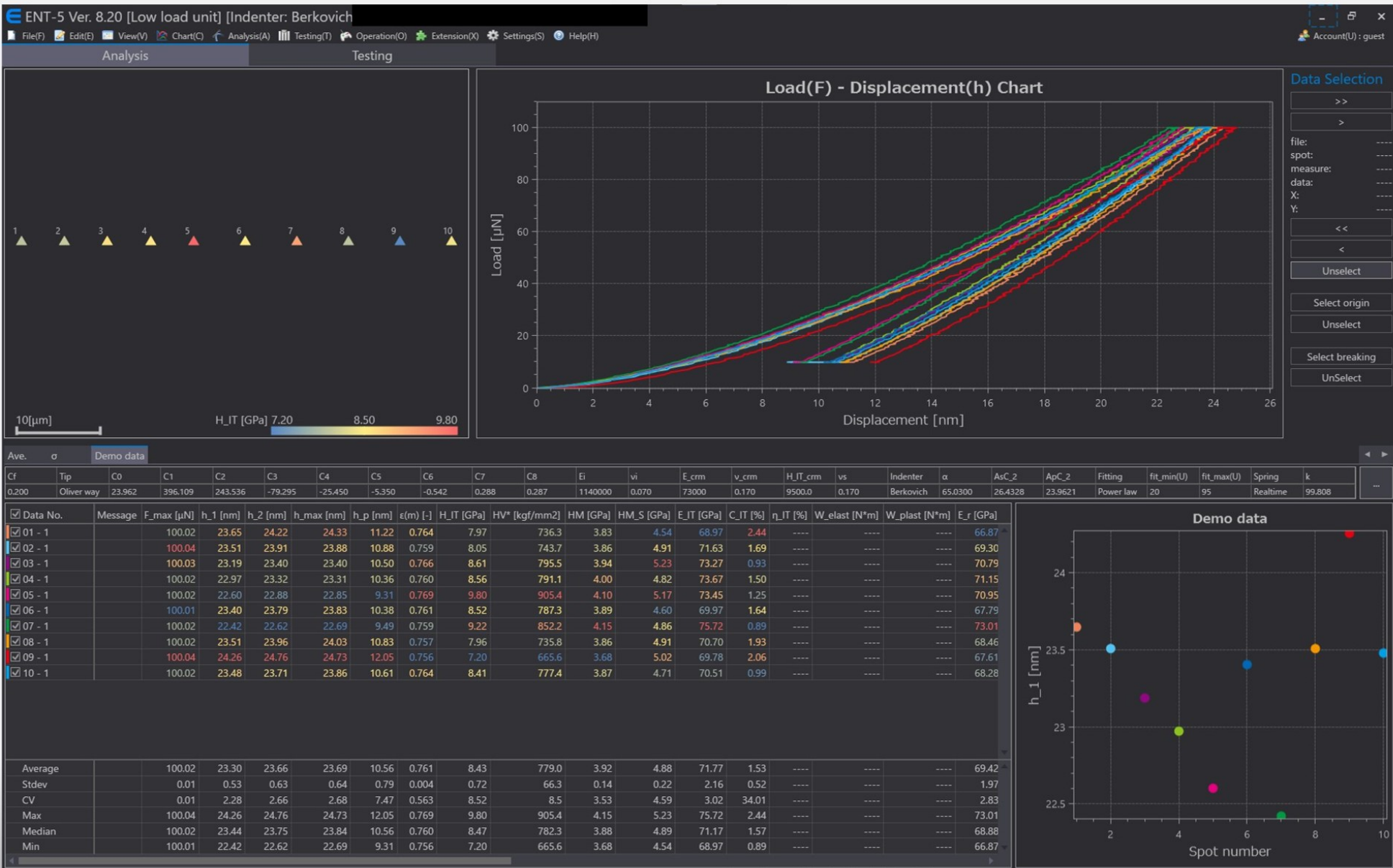

PDF report data can be
obtained by function of "PDF
out put".

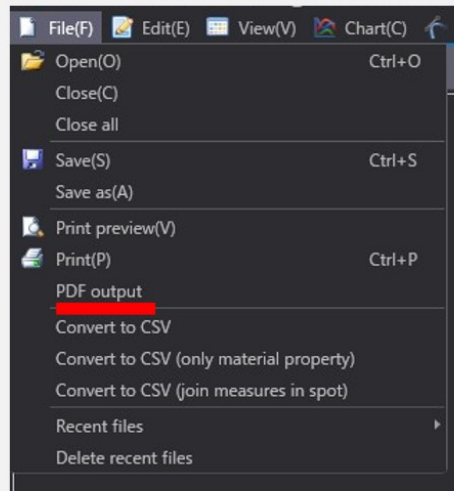

Video

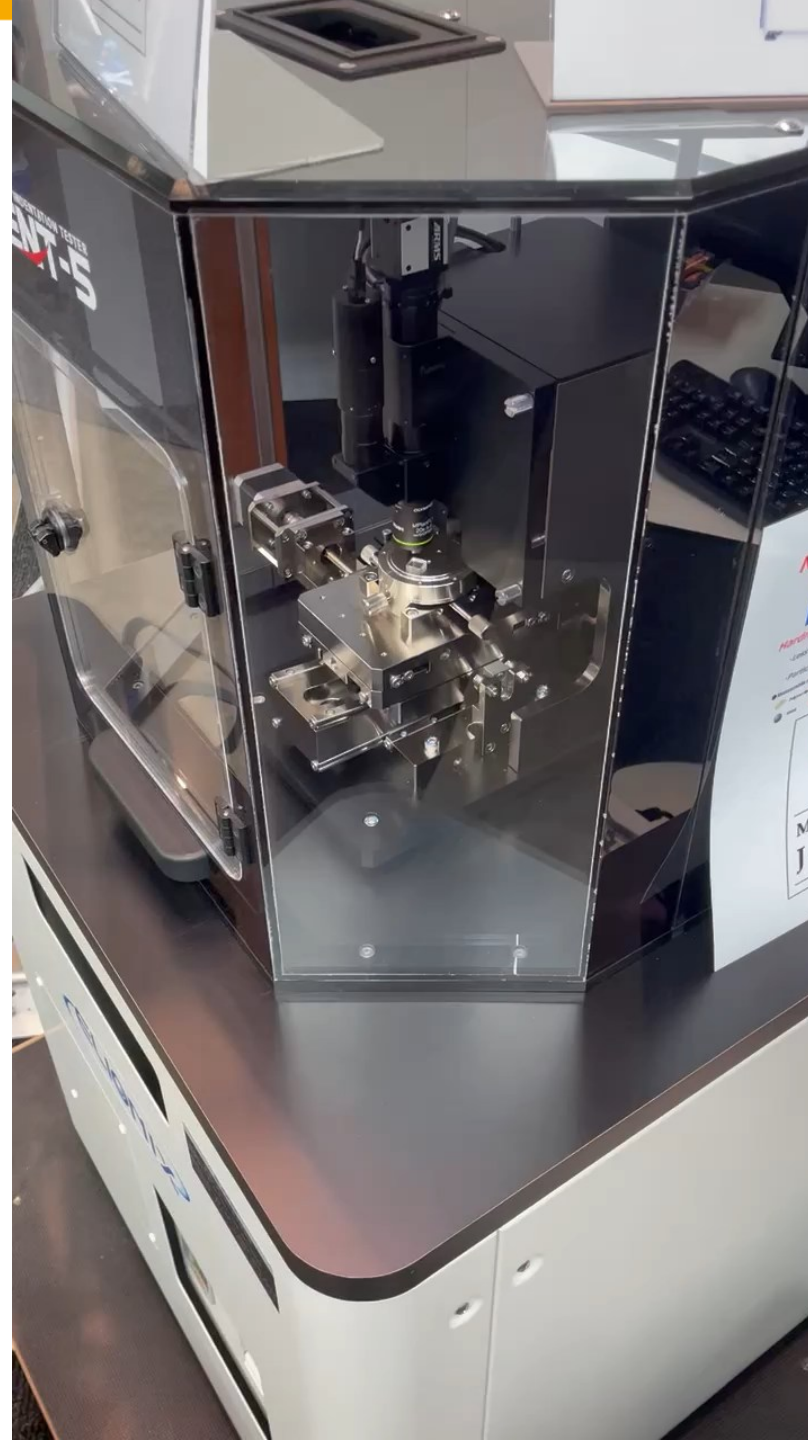



Dotazy

Děkuji za pozornost

Ing. Štěpán - Adam Havlíček

Product Specialist

Pragolab s.r.o.

havlicek@pragolab.cz

+420 771 274 373



Kontakt

